((SKY))

VIBRATION PLATE

SVP04

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EN MANUAL INSTRUCTION

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www.abisal.pl

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The ensure proper operation of the Product, please read this USER MANUAL carefully before using, and pay special attention to the safety precautions.

Contents

| | |
|-------------------------------------|----|
| Safety Precautions. | 8 |
| Function of the Product. | 9 |
| Programs. | 10 |
| Function of Controller. | 11 |
| Using Method. | 12 |
| Moving Method. | 13 |
| Cleaning and Maintenance. | 13 |
| Common Problems and Solutions. | 13 |
| Technical Specification. | 13 |

Safety Precautions

Environment for use

Do not use the Product in humid or dusty environment to avoid electric shock or breakdown of the mechanical parts.

Do not use the Product when indoor temperature is over 40°C.

Do not use the Product in the open air or under the direct sunlight.

Please the Product on flat floor.

Please wait for one hour before using the Product when it is moved from low to high temperature, so as to prevent water vapor from condensing on the machine and resulting in any failure.

People under these conditions below are not suggested to use the Product

It is strongly advised to consult your doctor before use

Patients who are under treatment or feel not good;

Patients with malignant tumor;

Pregnant women or women in menstrual period;

People with heart disease;

Osteoporosis patients;

People with heart pacemaker or other medical electronic instrument implanted;

People with artificial limb;

People who are sensitive to vibration and strenuous exercise.

Notices during Operation

Before use, please undress watch, bangle, call phone or other valuable objects, and keys, knife or other incisive objects, so as to avoid damage or body injury.

Do not smoke while using the Product.

Please do not use the Product together with other medical equipments.

Please do not use the Product for medical treatment purpose.

Do not use the Product right after a meal or drink.

During use, please keep child or pet away from the Product. It is forbidden to put fingers or other foreign bodies into the gap between the middle cover and base cover.

Supervision is needed for child under 10 years and the aged over 60 years during using the Product.

The load capability of the Product is 150 KGS; people over this weight are not suggested to use the Product.

This product is only limited for one- person use.

WARNING

Unauthorized dismantle, disassembly or renovation is not allowed to prevent electrical shock or any malfunctions. Please contact manufacture or local dealer immediately if there is any service work needed.

Stop using the Product under the Product under below conditions, then power off and remove the plug from power supply.

Water or other liquid is accidentally splashed onto the Product or controller.

Any accessories found damaged or internal component is exposed.

Feel abnormal or severe pain during use.

Power failure.

Thundering or lighting weather.

Plug and power cord

Check if the voltage complies with the one specified for the Product. Do not retract the power plug with wet hand to prevent electric shock. Hold the plug rather than the power cord to pull out the power cord. Turn off the power switch and remove the plug from the power supply in case of electrical failure or not in use for a long time. Stop using the Product if the power card has any damage, turn off the over switch, pull out the plug, and contact technicians for repair.

Function of the Product

P0 - abdomen stretching

P1 - relax stretching

P2 - gentle movement stretching

P3 - running stretching

P4 - jogging stretching

P5 - fast walking stretching

P6 - walking stretching

P7 - average stretching

P8 - high level stretching


P9 - heart stretching

Auto program speed sheet

Time from 15 min to 1 min

| | | 15min | 14min | 13min | 12min | 11min | 10min | 9min | 8min | 7min | 6min | 5min | 4min | 3min | 2min | 1min |
|------|----|-------|-------|-------|-------|-------|-------|------|------|------|------|------|------|------|------|------|
| low | PO | 4 | 8 | 10 | 12 | 17 | 20 | 24 | 20 | 17 | 6 | 12 | 17 | 20 | 24 | 20 |
| | P1 | 2 | 4 | 6 | 8 | 10 | 13 | 15 | 12 | 9 | 8 | 8 | 10 | 13 | 15 | 12 |
| | P2 | 8 | 14 | 8 | 14 | 8 | 20 | 8 | 14 | 8 | 14 | 14 | 8 | 20 | 8 | 14 |
| mid | P3 | 12 | 23 | 29 | 38 | 49 | 56 | 68 | 76 | 48 | 29 | 38 | 49 | 56 | 68 | 76 |
| | P4 | 9 | 13 | 16 | 21 | 24 | 35 | 67 | 50 | 38 | 30 | 21 | 24 | 35 | 67 | 50 |
| | P5 | 23 | 55 | 23 | 55 | 23 | 55 | 23 | 55 | 23 | 55 | 55 | 23 | 35 | 23 | 55 |
| | P6 | 11 | 11 | 25 | 25 | 38 | 38 | 48 | 48 | 25 | 25 | 25 | 38 | 38 | 48 | 48 |
| high | P7 | 30 | 30 | 50 | 50 | 70 | 70 | 40 | 40 | 60 | 60 | 50 | 70 | 70 | 40 | 40 |
| | P8 | 55 | 40 | 55 | 40 | 50 | 77 | 80 | 65 | 65 | 99 | 40 | 50 | 77 | 80 | 65 |
| | P9 | 40 | 85 | 75 | 65 | 25 | 96 | 78 | 60 | 65 | 80 | 65 | 25 | 96 | 78 | 60 |

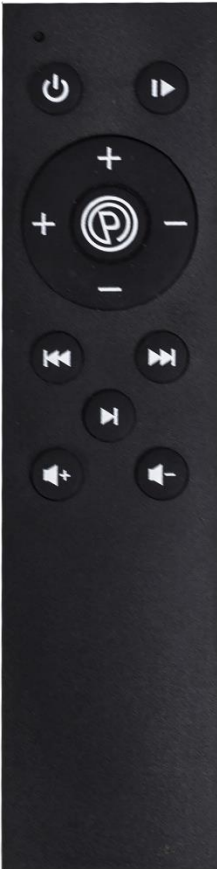START/ STOP: Insert the power plug, switch on power, then press the button to start the product. In the operation of the machine, press this button to stop working.

TIME-: Time decrease: the default working time is 10 minutes, divided into 10 levels, press this button to decrease time.

Speed+: Speed increase: Manual 1-20. Press this button to increase the speed when the machine is operating.

TIME+: Time increase: the default working time is 10 minutes, divided into 10 levels, press this button to increase time.

SPEED-: Speed decrease: Manual 1-20. Press this button to decrease the speed when the machine is operating.

Note: The effective reception range of the remote control signal is 2,5 meters, the infrared emission indicator of the controller should be aligned the infrared receiving window on the product.

1. Start: See 1
2. Program: See 2
3. Speed+/-: While machine works in manual mode, speed can be adjusted by this button.
4. Time+/-: When machine stop in manual mode, time can be adjusted between 1-20 minutes(Default time: 10minutes).
5. Led displaying method:
Method 1: Mode. Speed
Method 2: Minutes: Second
Both method will alternate in 5 seconds.

Use method

| | | | |
|---|---|---|--|
|  |  |  |  |
| <p>WHOLE-BODY-EXERCISE POSTURE: Stand on the plate and stretch your feet like as wide as your shoulders. This will help improve overall circulation and muscle tone.</p> | <p>SQUAT POSTURE: Stretch your legs as wide as your shoulders when squatting and bend your knees down to 90 degrees if you can and hold the position. This is great for a thigh workout.</p> | <p>PRESS-UP POSTURE: Stretch your arms as wide as your shoulders. If you wish to exercise the muscles more, bend your elbows to 90 degrees and hold this position.</p> | <p>LOWER LEGS POSTURE: Put your lower legs as shown on the plate with your hands supporting behind you on the ground. For a deeper exercise, raise your seat off the floor and press down on your calves.</p> |
|  |  |  |  |
| <p>WAIST BEND POSTURE: Stretch your legs, bend at the waist and put your hands on the plate. Bend your elbows for a deeper exercise.</p> | <p>SINGLE-FOOT POSTURE: Put one foot on the pedal and relax yourself. For a deeper exercise, lunge forward, placing more body weight on the plate.</p> | <p>SITTING POSTURE I: Sit with your buttocks on the oscillating plate. Sit comfortably, keeping your back straight.</p> | <p>SITTING POSTURE II: Put your legs on the plate when sitting on a chair. This is a great way to ease tension and can help to improve circulation in the legs.</p> |


Function and features

Body beauty: by reposeful and high-frequency vibration and sway, combined with various exercise postures, this machine can deeply stimulate muscle tissue, lower fat accumulation, rehabilitate elasticity of muscle, sculpture beautiful body shape, then make users vitalized and feel confident.

Exercise: Borrow health and fitness ideas from Chinese traditional culture like Kung Fu, Buddhism, Taoism, it advocates passive exercise, which allow you to exercise indoor in a relatively static manner to get mind and body relaxed, pain and soreness relived, fatigue swept away, inner circulation improved, then make you stronger and healthier.

Moving method

Handle - Two handles are designed on both ends of the middle cover. When moving the Product, lift the end opposite of the castor, grasp the handle and move the Product to the desired place.

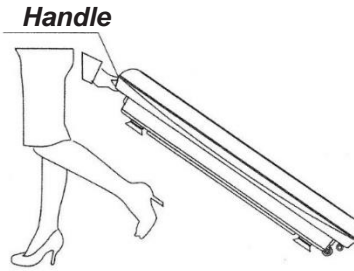Storage and maintenance surface cleaning

Before cleaning, ensure the power switch is turned off, and the attaching plug is pulled out of the socket. The rubber pedal and covers can be cleaned with a wet cloth containing detergent, and the detergent should be removed with a cloth.

NOTICE: Do not clean the product with corrosive detergent such as benzene or thinner, as it may cause fading, corrosion or crack on the Product. Do not splash water into the product to avoid short circuit.

CONTROLLER

As wet cloth may cause malfunction of the controller, only dry cloth can be used.

Storage and maintenance

Make sure that the power switch is turned off and the plug is removed from the power supply. Clean off the dust and dirt on the Product and keep it clean and away from places with high temperature or high humidity.

If the Product is not in use for a long period, a dustproof cover (but it yourself) should be applied on the Product to prevent the entry of dust.

Common problems & Solutions

1. Emergency stop: Led on console shows "Er-3", Led on PCB flashes three Times.
2. Connection error: Led on console shows "ER-4", Led on PCB flashes four Times.
3. Electron power tube damage: Led on console shows "ER-5", Led on PCB flashes five Times.
4. Current Overload: Led on console shows "ER-6", Led on PCB flashes six Times.
5. Engineer abnormality: Led on console shows "ER-7", Led on PCB flashes seven Times.
6. Blocked Engineer: Led on console shows "ER-8", Led on PCB flashes eight Times.
7. Low voltage: Led on console shows "ER-9", Led on PCB flashes nine Times.

Technical specification

Rated voltage: 220-240V

Rated power: 240W Rated frequency: 50/60 Hz Maximum user weight: 150 kg

Assembly size: 69*39*13cm

G.W./N.W: 14,5/ 12,5kg, **Functional massage with magnetic balls**

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Article name: EAN code: Date of sale:

GUARANTEE TERMS

- The Seller provides guarantee on behalf of the Guarantor within the territory of the Republic of Poland for the period of 24 months from the date of sale.
 - The Guarantee will be recognised by the shop or service centre after the client provides:
 - clearly and correctly filled-in guarantee card with the sale stamp and the seller's signature
 - valid purchase confirmation for the equipment including the date of sale / bill /
 - complained product.
 - Possible defects and damages revealed during the guarantee period shall be remedied free of charge within not more than 21 days from the delivery of the product to the shop or service.
 - Should it be necessary to import some components, the guarantee term can be extended by the period required to import such part, however such period shall not be longer than 40 days.
 - The guarantee does not cover:
 - mechanical damages and subsequent defects,
 - damages and defects resulting from incorrect storage and usage of the equipment against its purpose,
 - incorrect assembly and maintenance,
 - damages and consumption of such elements as: cables, straps, rubber parts, pedals, sponge grips, wheels, bearings, etc.
 - The guarantee expires in the event of:
 - expiry of the validity term,
 - self-repairs,
 - failure to follow the rules of correct operation.
 - No duplicates of the guarantee cards shall be issued.
 - The product handed over for repair should be complete and clean. If some parts are missing, the service is entitled to refuse to accept the product for repair. The service may refuse to accept a dirty product or clean in on the client's expense after the client's written permit.
 - The guarantee does not cover activities related to assembly, maintenance, which, according to the instruction, the user is obliged to perform themselves.
 - The Guarantor also informs that they offer post-guarantee service.
 - The product should be provided in its original packaging and should be secured for shipment.
- The guarantee for sold goods is not exclusive, does not limit or suspend entitlements of the buyer resulting from the non-compliance of the product with the agreement.

THE EQUIPMENT IS NOT INTENDED FOR REHABILITATION AND THERAPY

NOTES ON THE COURSE OF REPAIRS

| Item | Date of notification | Date of provision | Course of repairs | Signature of the recipient (shop, owner) |
|------|----------------------|-------------------|-------------------|--|
| | | | | |
| | | | | |
| | | | | |
| | | | | |

Even though we do our best to ensure the highest quality of our products, single defects or omissions might occur. If you note a defect or missing part, we invite you to contact us.

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Labelling the product with a symbol of a crossed waste container informs us that the worn-out electrical and electronic equipment cannot be put into the same container as other waste. According to the WEEE Directive on the management of electrical and electronic waste, separate recycling methods shall be used for this type of equipment.

A user who wants to dispose of such product is obliged to transfer it to the collection centre for worn-out electrical and electronic equipment, thanks to which they support its re-usage, recycling or recovery and thusly protect natural environment. For this purpose, one should contact the equipment seller or representatives of local authorities. Hazardous components contained in electronic equipment may have long-lasting harmful impact on natural environment as well as adverse effect on human health.

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